

HELIOS 785

BEFORE AND AFTER PHOTO

LASEROPTEK
Amico Yasna Pars (Pr.J.S.Co)

MELASMA LIGHT PIGMENTS FULL FACE TONING

1064 Zoomcollimator, 8mm, 1.4J
1064 Zoomcollimator, 7mm, 2.3J
785 FX, 5x5mm, 0.8J
785 Zoom, 4mm, 1.0J, 3mm, 2.0J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



BEFORE

MELASMA LIGHT PIGMENTS FULL FACE TONING

1064 Zoomcollimator, 8mm, 1.4J
1064 Zoomcollimator, 7mm, 2.3J
785 FX, 5x5mm, 0.8J
785 Zoom, 4mm, 1.0J, 3mm, 2.0J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE

AFTER 5TX

AFTER 10TX

MELASMA LIGHT PIGMENTS FULL FACE TONING

1064 Zoomcollimator, 8mm, 1.4J
1064 Zoomcollimator, 7mm, 2.3J
785 FX, 5x5mm, 0.8J
785 Zoom, 4mm, 1.0J, 3mm, 2.0J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE

AFTER 5TX

AFTER 10TX

MELASMA LIGHT PIGMENTS FULL FACE TONING

1064 Zoomcollimator, 8mm, 1.4J
1064 Zoomcollimator, 7mm, 2.3J
785 FX, 5x5mm, 0.8J
785 Zoom, 4mm, 1.0J, 3mm, 2.0J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

MELASMA LIGHT PIGMENTS FULL FACE TONING

1064 Zoomcollimator, 8mm, 1.4J
1064 Zoomcollimator, 7mm, 2.3J
785 FX, 5x5mm, 0.8J
785 Zoom, 4x4mm, 1.0J, 3mm, 2.0J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

MELASMA

785 FX, 5x5mm, 7-10hz, 0.4J, 1-2passes
1064 Zoomcollimator, FR, 8mm, 10hz, 4.0J, 1-2 passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO PIH REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes
785 Zoom, 3x3mm, 2.0J
1064 Zoomcollimator, FR, 8mm, 10hz, 5.0J-J, 1-2 passes,

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

SKIN WHITENING BLEMISHES REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.8J, 1-2passes
1064 Zoomcollimator, FR, 8mm, 5.0-5.5J, 10hz, 1-2 passes
1064 Zoomcollimator, 8mm, 1.3J, 1-2passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

SKIN WHITENING LENTIGO REJUVENATION SKIN BOOSTER

785 FX, 5x5mm, 7-10hz, 0.6-0.8J, 2-3passes
Combination therapy – Skin booster

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

SKIN WHITENING LENTIGO REJUVENATION SKIN BOOSTER

785 FX, 5x5mm, 7-10hz, 0.6-0.8J, 2-3passes
Combination therapy – Skin booster

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

SKIN WHITENING LENTIGO REJUVENATION SKIN BOOSTER

785 FX, 5x5mm, 7-10hz, 0.6-0.8J, 2-3passes
Combination therapy – Skin booster

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

SKIN WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.7-0.8J, 2-3passes
1064 Zoomcollimator, 5mm, 4.0-4.6J, 2-3passes,
6TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

SKIN WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.7-0.8J, 2-3passes
1064 Zoomcollimator, 5mm, 4.0-4.6J, 2-3passes,
6TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

SKIN WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.7-0.8J, 2-3passes
1064 Zoomcollimator, 5mm, 4.0-4.6J, 2-3passes,
6TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE




AFTER

SKIN WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.7-0.8J, 2-3passes
1064 Zoomcollimator, 5mm, 4.0-4.6J, 2-3passes,
6TX

Amico Yasna Pars (Pr.J.S.Co)

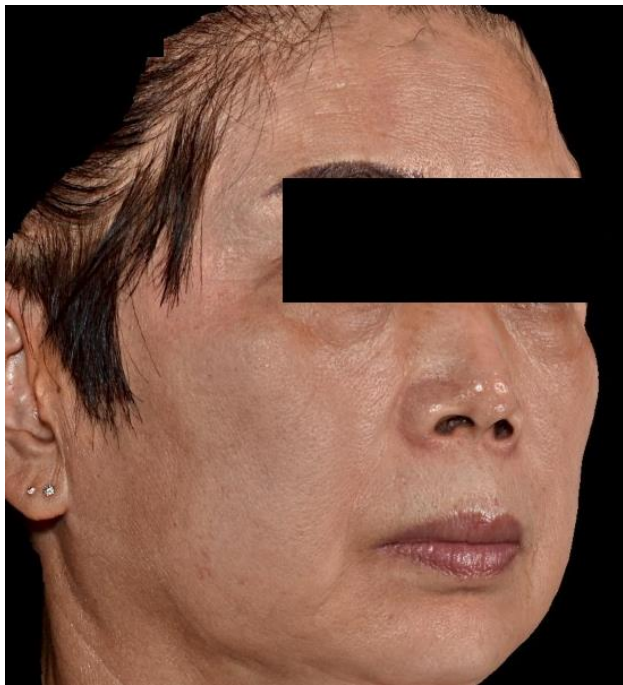

BEFORE



AFTER


SKIN WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.7-0.8J, 2-3passes
1064 Zoomcollimator, 5mm, 4.0-4.6J, 2-3passes,
5TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE

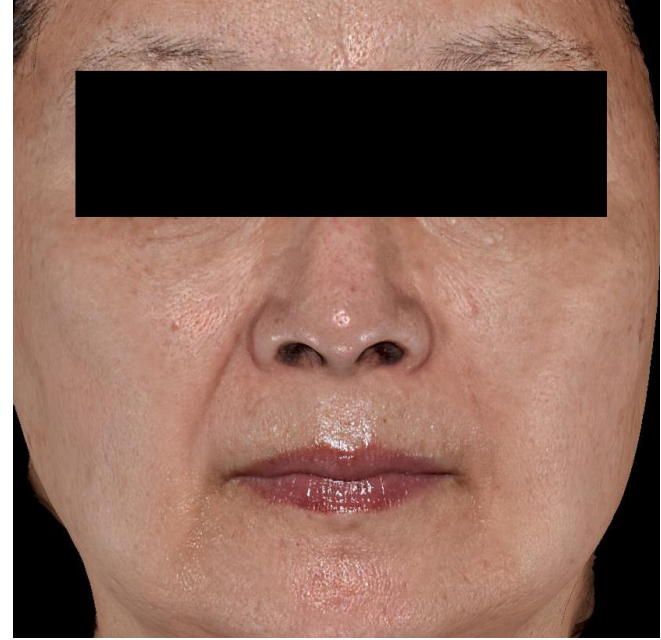

AFTER

SKIN WHITENING LENTIGO SKIN TEXTURE

785 FX, 5x5mm, 7-10hz, 0.7-0.8J, 2-3passes
785 Zoom, 3x3mm, 1-3hz, 1.2-1.5J
Skin texture – combination therapy using Lotus III

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO SKIN WHITENING FULL FACE TONING

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes
785 Zoom, 4x4mm, 1.2J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO SKIN WHITENING FULL FACE TONING

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes
785 Zoom, 4x4mm, 1.2J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER


LENTIGO SKIN WHITENING FULL FACE TONING

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes
785 Zoom, 4x4mm, 1.2J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE

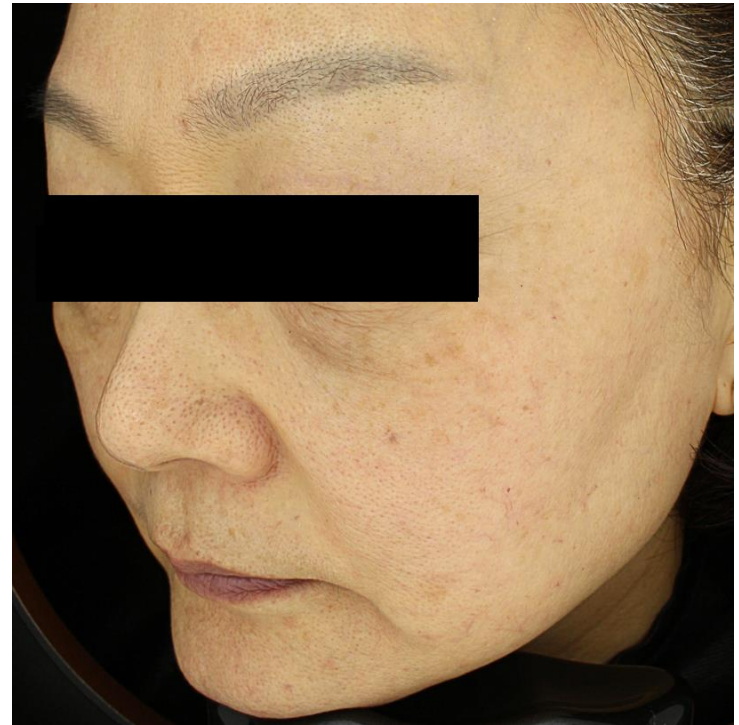


AFTER

LENTIGO WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes
785 Zoom, 3x3-4x4mm, 1.2-1.8J

Amico Yasna Pars (Pr.J.S.Co)

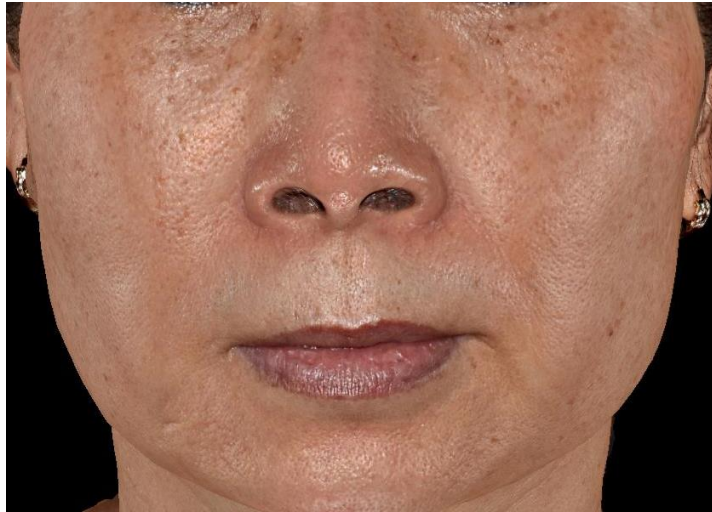

BEFORE



AFTER

LENTIGO WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes
785 Zoom, 3x3-4x4mm, 1.2-1.8J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO WHITENING REJUVENATION

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes
785 Zoom, 3x3-4x4mm, 1.2-1.8J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO

785 Zoom, 3x3-4x4mm, 1.2-1.5J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO

785 Zoom, 3x3-4x4mm, 1.2-1.5J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO

785 Zoom, 3x3-4x4mm, 1.2-1.5J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO SKIN TEXTURE

785 Zoom, 3x3-4x4mm, 1.2-1.5J
785 FX, 5x5, 7-10hz, 0.8J, 2-3passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

FRECKLE SKIN TEXTURE

785 Zoom, 3x3-4x4mm, 1.2-1.5J
785 FX, 5x5, 7-10hz, 0.7-0.8J, 2-3passes, 7TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

FRECKLE

785 Zoom, 4x4mm, 1.2J, 1TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO 785 TONING

785 Zoom, 4x4mm, 1.0-1.2J
785 FX, 5x5mm, 0.8J, 2-3passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO 785 TONING

785 Zoom, 4x4mm, 1.0-1.2J
785 FX, 5x5mm, 0.8J, 2-3passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO 785 TONING

785 Zoom, 4x4mm, 1.0-1.2J
785 FX, 5x5mm, 0.8J, 2-3passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE




AFTER

LENTIGO

785 Zoom, 4x4mm, 1.0-1.2J, 1TX

Amico Yasna Pars (Pr.J.S.Co)

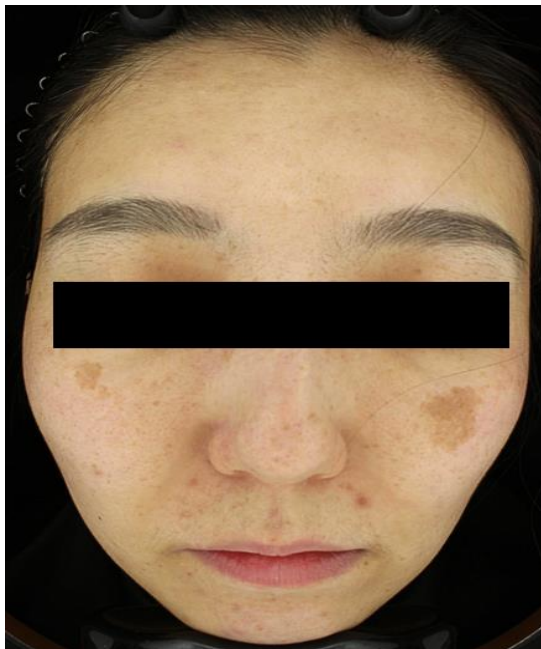

BEFORE



AFTER

FRECKLE

785 Zoom, 3X3mm, 1.5J, 1TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

LENTIGO SKIN WHITENING

785 Zoom, 3x3mm, 1.3J
785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE




AFTER

CENTROFACIAL LENTIGINOSIS

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes, 5TX

Amico Yasna Pars (Pr.J.S.Co)

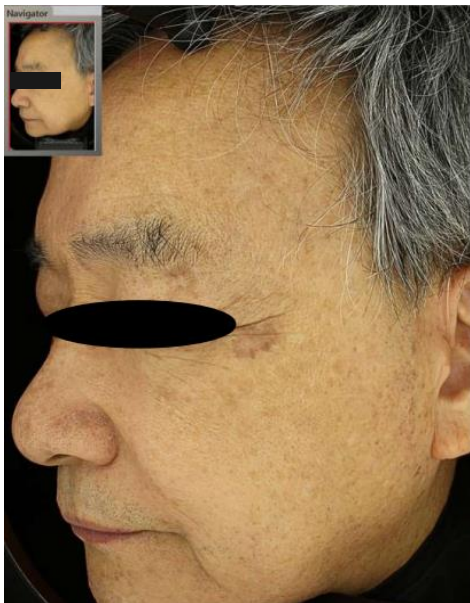

BEFORE



AFTER 2TX



AFTER 5TX

FRECKLE LENTIGO

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes, 9TX
785 Zoom, 3x3mm, 1.4J, 2-5hz

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

FRECKLE LENTIGO

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes, 9TX
785 Zoom, 3x3mm, 1.4J, 2-5hz

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

FRECKLE LENTIGO

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes, 9TX
785 Zoom, 3x3mm, 1.4J, 2-5hz

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

ABNOM LENTIGO

785 Zoom, 4x4mm, 1.2J
785 FX, 5x5mm, 0.8J, 2-3passes, 8TX
A case after 1064 toning

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



**BEFORE
(Q-1064 Toning 10TX)**



AFTER (P-785 8TX)

ACNE PIH WHITENING REJUVENATION SKIN TEXTURE

785 FX, 5x5mm, 0.8J, 2-3passes
785 Zoom, 4x4mm, 1.2J

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

ERYTHROMELANOSIS FOLLICULARIS FACIAL ET COLLI

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

RIEHL'S MELANOSIS

785 FX, 5x5mm, 7-10hz, 0.8J, 2-3passes, 7TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

BECKER'S NEVUS

785 FX, 5x5mm, 10hz, 0.8J, 5000 shots
1064 Zoomcollimator, 8mm, 10hz, 2.7J, 9TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER

TATTOO REMOVAL

785 Zoom, 3x3-5x5mm, 0.8-1.9J, 2hz, 1pass
1064 Zoom 4-5mm, 5.0-7.0J, 2hz, 1pass
532 Zoom 3-4mm, 2.0-2.7J, 2hz, 1pass, 4TX

Amico Yasna Pars (Pr.J.S.Co)



BEFORE



AFTER