

LASEROPTTEK

Amico Yasna Pars (Pr.J.S.Co)

PICOLO PREMIUM

BEFORE AND AFTER PICTURES BOOK

Scars treatment

- **Automutilation scars**
- 1st Fractional 1064 HP I (10x10mm), 0.30J/cm², II 0.20J/cm², III 0.13J/cm², 7-10hz
- 2ND Fractional 1064 HP I (10x10mm), 0.40J/cm², II 0.30J/cm², III 0.15J/cm², 7-10hz
- 3rd Fractional 1064 HP I (10x10mm), 0.40J/cm², II 0.30J/cm², III 0.14J/cm², 7-10hz



Before



After 3Tx

- Acne scar
- Fractional 1064 HP (10x10mm), 0.26-0.30J/cm², 5-10hz, 3-5passes



Before



After 3 Tx

- Lower Dermal layer : Fractional 1064 HP I (10x10mm), 0.5J/cm², 5-10hz, 3passes 5000 shots
- Upper dermal layer : Fractional 1064 HP II (10x10mm), 0.4J/ cm², 5-10hz, 2passes 3000 shots
- Epidermal layer : Fractional 1064 HP III (10x10mm), 0.36J/cm², 5hz, 3passes 2000 shots



Before



After 6 Tx

- Lower Dermal layer : Fractional 1064 HP I(10x10mm), 0.5J/cm², 5-10hz, 3passes 5000 shots
- Upper dermal layer : Fractional 1064 HP II (10x10mm), 0.4J/ cm², 5-10hz, 2passes 3000 shots
- Epidermal layer : Fractional 1064 HP III (10x10mm), 0.36J/cm², 5hz, 3passes 2000 shots



Before



After 6 Tx

- Lower Dermal layer : Fractional 1064 HP I(10x10mm), 0.5J/cm², 5-10hz, 3passes 5000 shots
- Upper dermal layer : Fractional 1064 HP II (10x10mm), 0.4J/ cm², 5-10hz, 2passes 3000 shots
- Epidermal layer : Fractional 1064 HP III (10x10mm), 0.36J/cm², 5hz, 3passes 2000 shots



Before



After 6 Tx

- Lower Dermal layer : Fractional 1064 HP I (10x10mm), 0.5J/cm², 5-10hz, 3passes 5000 shots
- Upper dermal layer : Fractional 1064 HP II (10x10mm), 0.4J/ cm², 5-10hz, 2passes 3000 shots
- Epidermal layer : Fractional 1064 HP III (10x10mm), 0.36J/cm², 5hz, 3passes 2000 shots



Baseline

D1

D2

D3

- Lower Dermal layer : Fractional 1064 HP I(10x10mm), 0.5J/cm², 5-10hz, 3passes 5000 shots
- Upper dermal layer : Fractional 1064 HP II (10x10mm), 0.4J/ cm², 5-10hz, 2passes 3000 shots
- Epidermal layer : Fractional 1064 HP III (10x10mm), 0.36J/cm², 5hz, 3passes 2000 shots



Before



After 4 Tx

- Acne scar
- Fractional 1064 HP I (10x10mm), 0.3-0.5J/cm², 8-10hz, 3-5passes



Before



After 3Tx

- Acne scar
- Fractional 1064 HP I (10x10mm), 0.3-0.5J/cm², 8-10hz, 3-5passes



Before



After 3Tx

- Acne scar
- Fractional 1064 HP (10x10mm), 0.3-0.5J/cm², 8-10hz, 3-5passes



Before



After 3Tx

- Acne scar
- Fractional 1064 HP I (10x10mm), 0.3-0.5J/cm², 8-10hz, 3-5passes



Before



After 3Tx

- Scar with PIH
- Fractional 1064 HP (10x10mm), 0.20-0.28J/cm², 5-10hz, 3-5passes



Before



After

- Depressed scar
- Fractional 1064 HP III (10x10mm), 0.26-0.30J/cm², 5hz, 3-5passes
- Combination Therapy with needle RF



Before




After 4Tx

- Acne scar
- 1064 nm Dia FX HP, 4 sessions
- Level I _ $0.45\text{J}/\text{cm}^2$, 8-10hz
- Level II _ $0.3\text{J}/\text{cm}^2$, 7hz
- Level III _ $0.2\text{J}/\text{cm}^2$, 3hz



Before

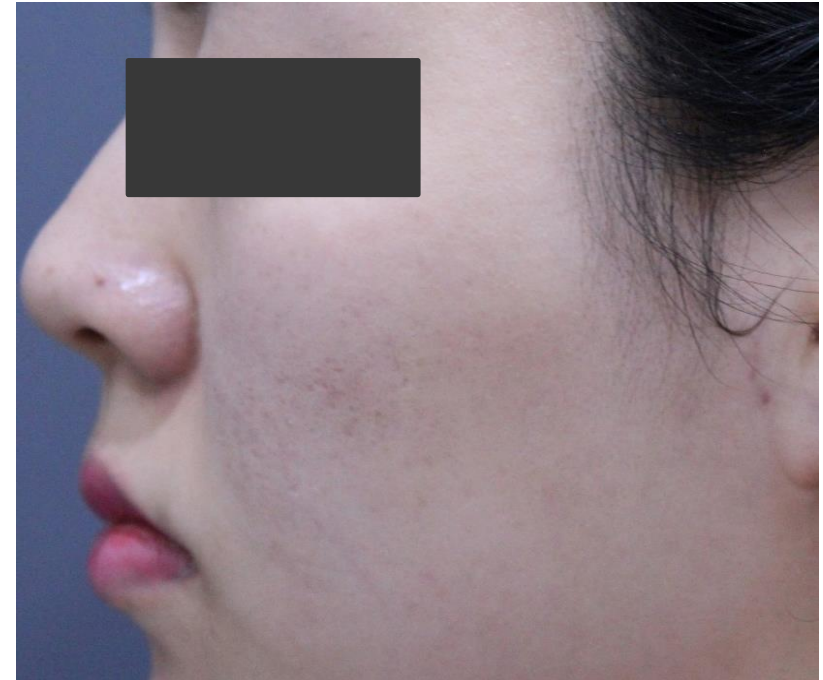

After 4 TX

- Acne scar
- 1064 nm Dia FX HP, 4 sessions
- Level I _ $0.5\text{J}/\text{cm}^2$, 8-10hz
- Level II _ $0.3\text{J}/\text{cm}^2$, 7hz
- Level III _ $0.2\text{J}/\text{cm}^2$, 3hz



Before



After 4 TX

- Acne scar
- 1064 nm Dia FX HP, 4 sessions
- Level I _ $0.45\text{J}/\text{cm}^2$, 8-10hz
- Level II _ $0.3\text{J}/\text{cm}^2$, 7hz
- Level III _ $0.2\text{J}/\text{cm}^2$, 3hz



Before



After 1 TX



After 2 TX



After 3 TX



After 4 TX

- Acne scar
- 1064 nm Dia FX HP, 4 sessions
- Level I _ $0.5\text{J}/\text{cm}^2$, 8-10hz
- Level II _ $0.3\text{J}/\text{cm}^2$, 7hz
- Level III _ $0.2\text{J}/\text{cm}^2$, 3hz



Before



After 4 TX

- Skin texture (rejuvenation, scars)
- 1064 nm Dia FX HP, 3 sessions per 8 weeks
- Level II, 10x10 mm, 3 hz, 0.24 J/cm²



- Acne scar
- 1064 nm Dia FX HP, 4 sessions
- Level II _ 0.38J/cm², 4hz



Before



After 4 TX

- Acne scar
- Fractional 1064 HP I (10x10mm), 0.3-0.5J/cm², 5-10hz, 3-5passes
- Fractional 1064 HP II (10x10mm), 0.15-0.20J/cm², 5hz, 3-5passes



Before



After

- Depressed scars, Enlarged pores
- Fractional 1064 HP I (10x10mm), 0.3-0.5J/cm², 5-10hz, 3-5passes
- Fractional 1064 HP II (10x10mm), 0.15-0.20J/cm², 5hz, 3-5passes



Before



After

- Enlarged pores
- Fractional 1064 HP (10x10mm), 0.38J/cm², 5hz, 5passes



Before



After 3Tx

Tattoo removal treatment

- 1st : Zoom HP, 1064nm, 4mm, 1.5J/cm², 1-3hz, 1pass
Fractional HP, 1064nm, Level I, 0.30J/cm², 5-10hz, 3passes
- 2nd : Zoom HP, 1064nm, 5mm, 1.5J/cm², 1-3hz, 1pass
Fractional HP, 1064nm, Level I, 0.40J/cm², 5-10hz, 4passes



Before



After 2Tx

- 1st : Zoom HP, 1064nm, 4mm, 1.5J/cm², 1-3hz, 1pass
Fractional HP, 1064nm, Level I, 0.30J/cm², 5-10hz, 3passes
- 2nd : Zoom HP, 1064nm, 5mm, 1.5J/cm², 1-3hz, 1pass
Fractional HP, 1064nm, Level I, 0.40J/cm², 5-10hz, 4passes




Before



After 2Tx

- 1st : QSNL
- 2nd : PIP 1064nm, zoom HP, 5mm, 1.7J/cm², 1hz, 1pass
- 3rd : PIP 1064nm, zoom HP, 6mm, 1.1J/cm², 1hz, 1pass

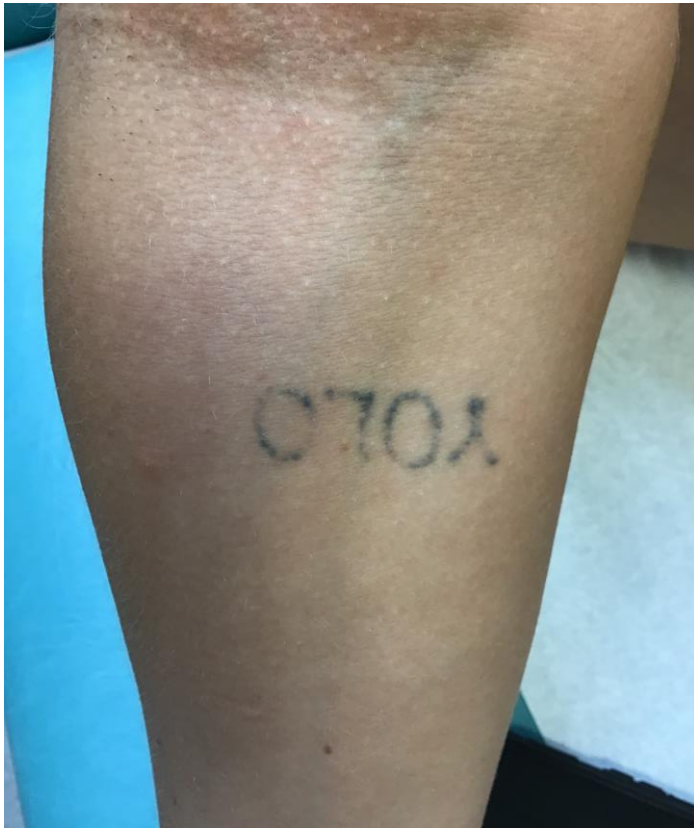

Before




After 2Tx

- Zoom HP, 1064nm, 5mm, 1.3-1.8J/cm², 1-3hz, 1pass



Before

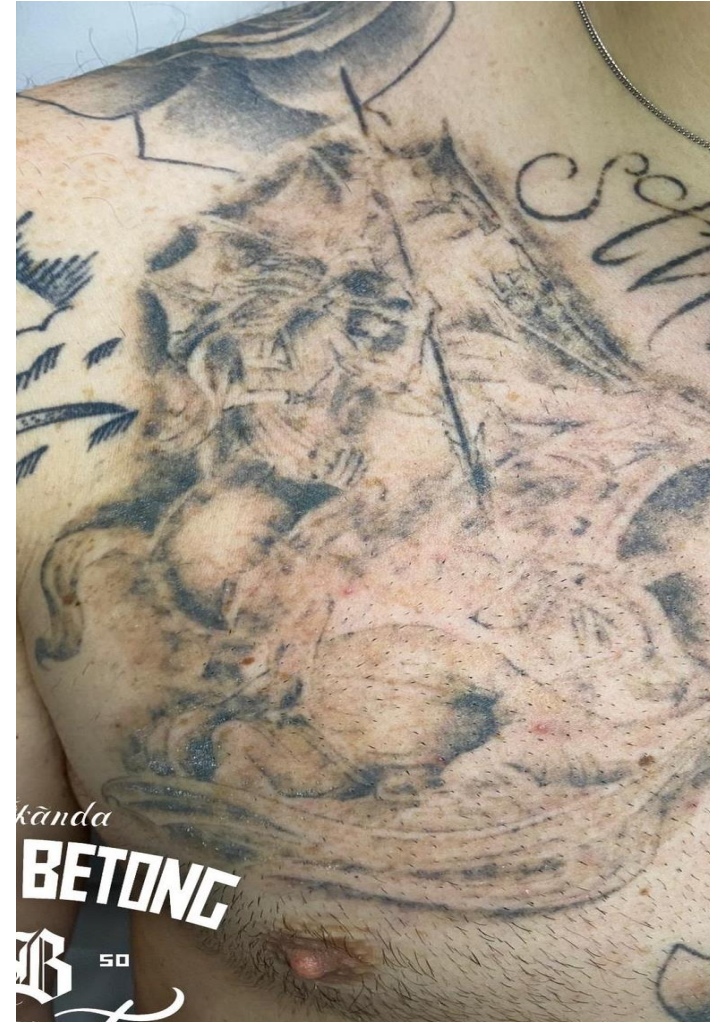

After 1Tx

- Zoom HP, 1064nm, 5mm, 1.3-1.8J/cm², 1-3hz, 1pass

Before



After 1Tx



- Zoom HP, 1064nm, 5mm, 1.3-1.8J/cm², 1-3hz, 1pass




Before



After 1Tx

- Zoom HP, 1064nm, 5mm, 1.3-1.8J/cm², 1-3hz, 1pass

Before

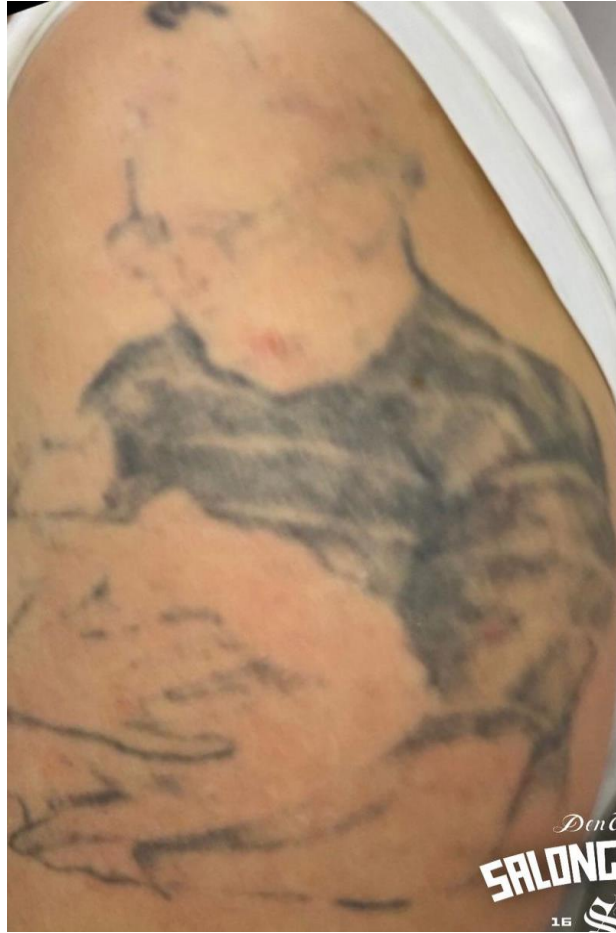

After 3 Tx



- Black : Zoom HP, 1064nm, 5mm, 1.3-1.8J/cm², 1-3hz, 1pass
- Red, orange n yellow : Zoom HP, 532nm, 5mm, 0.6-0.9J/cm², 1-3hz, 1pass



Before



After 1Tx

- Zoom HP, 1064nm, 3mm, 1.9J/cm², 1-3hz, 1pass



Before




After 5Tx

- Zoom HP, 3mm, 1064nm, 2-3J/cm², 1-3hz, 1pass



Before

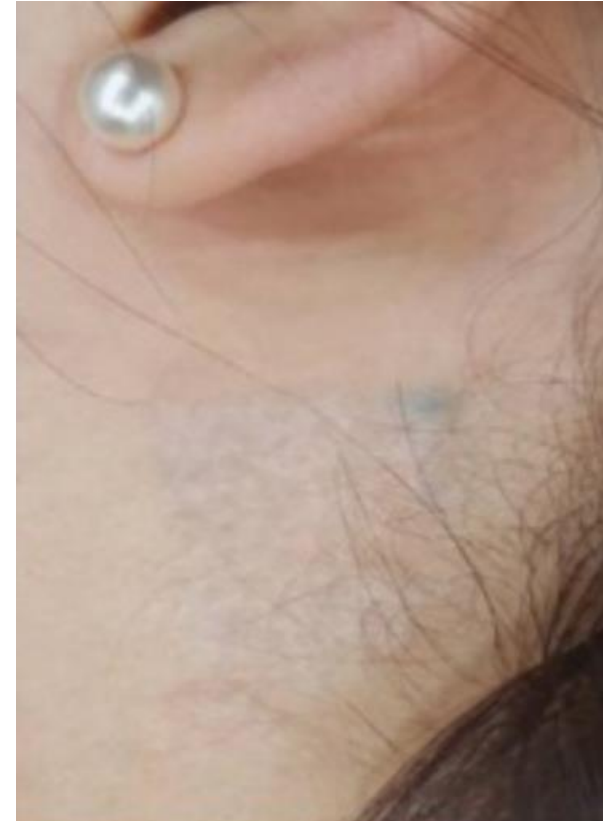

After 5Tx

- Black : Zoom HP, 4-5mm, 1064nm, 1.3-1.6J/cm², 1-3hz, 1pass
- Yellow, Red : Zoom HP, 4-5mm, 532nm, 0.4-0.6J/cm², 1-3hz, 1pass



Before

After 3Tx

- Black colored tattoo removal
- 1064 nm Zoom HP, 3 sessions
- 5mm 1.5-1.8 J/cm²



Before



After 3 TX

- Black colored tattoo removal
- 1064 nm Zoom HP, 3 sessions
- 5mm 1.5-2.0 J/cm²




Before



After 3 TX

- Black colored tattoo removal
- 1064 nm Zoom HP, 5 sessions
- 5mm 3-5HZ, 1.5-1.8 J/cm² → 4mm 1.8-2.0 J/cm²

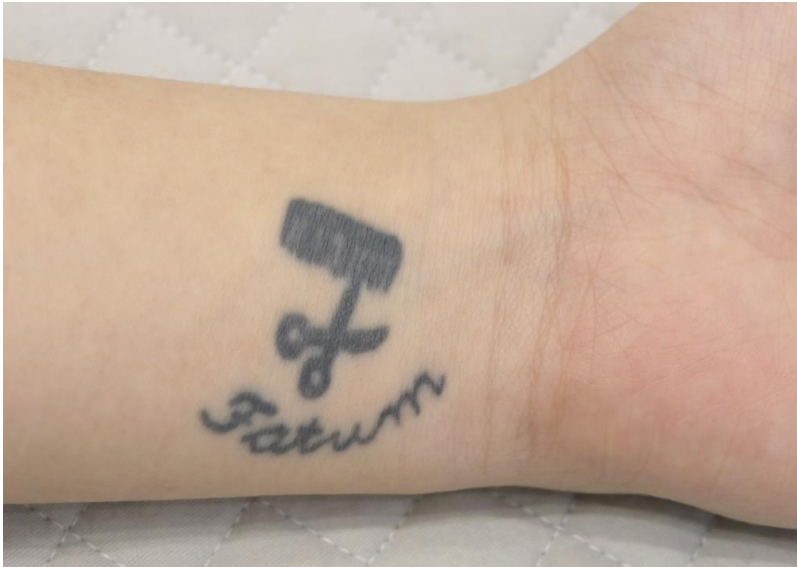

Before



After 5 TX

- Black colored tattoo removal
- 1064 nm Zoom HP, 4 sessions
- 5mm 3-5 HZ, 1.5-1.8 J/cm²



Before




After 4 TX

- Black colored tattoo removal
- 1064 nm Zoom HP, 4 sessions
- 3-5 HZ, 5mm 1.7- 2.2 J/cm² → 4mm 1.8-2.0 J/cm²



Before

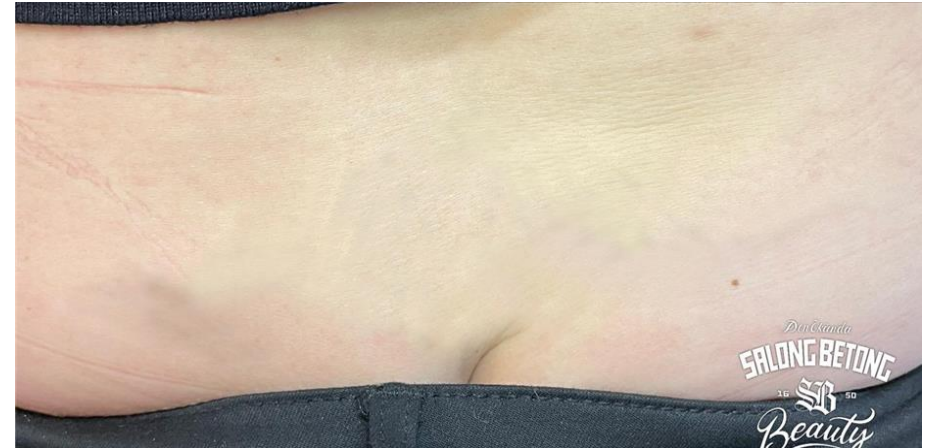

After 5 TX

Pico toning, whitening and skin rejuvenation treatment

- Jawline Contour, skin tightening, skin rejuvenation
- Tx1 _ 1064 Dia FX I (10x10mm), 0.2J/cm², 7hz, 8,000 shots to the face to under chin, multiple passes
- Tx2 _ 1064 Dia FX I (10x10mm), 0.3J/cm², 7hz, 10,000 shots, face to under chin, moving lower onto neck
- Tx3 _ 1064 Dia FX I (10x10mm), 0.4J/cm², 8hz, 10,000 shots, face to under chin, moving lower onto neck
- Tx4 _ 1064 Dia FX I (10x10mm), 0.4J/cm², 6hz, 12,000 shots for full face, 4000 shots for pore



- Jawline Contour, skin tightening, skin rejuvenation
- Tx1 _ 1064 Dia FX I (10x10mm), 0.2J/cm², 7hz, 8,000 shots to the face to under chin, multiple passes
- Tx2 _ 1064 Dia FX I (10x10mm), 0.3J/cm², 7hz, 10,000 shots, face to under chin, moving lower onto neck
- Tx3 _ 1064 Dia FX I (10x10mm), 0.4J/cm², 8hz, 10,000 shots, face to under chin, moving lower onto neck
- Tx4 _ 1064 Dia FX I (10x10mm), 0.4J/cm², 6hz, 12,000 shots for full face, 4000 shots for pore



- Skin Rejuvenation, Skin Tightening, Skin Texture
- Fractional 1064 HP I (10x10mm), 0.30-0.38J/cm², 5-10hz, 3-5passes
- Fractional 1064 HP II (10x10mm), 0.20-0.25J/cm², 5-10hz, 3-5passes



Before



After

- Skin Rejuvenation, Skin Tightening, Skin Texture
- Fractional 1064 HP I (10x10mm), 0.30-0.38J/cm², 5-10hz, 3-5passes
- Fractional 1064 HP II (10x10mm), 0.20-0.25J/cm², 5-10hz, 3-5passes



Before



After

- Skin Rejuvenation, Skin Tightening, Skin Texture
- Fractional 1064 HP I (10x10mm), 0.30-0.38J/cm², 5-10hz, 3-5passes
- Fractional 1064 HP II (10x10mm), 0.20-0.25J/cm², 5-10hz, 3-5passes



Before



After

- Skin Rejuvenation, Skin Tightening, Skin Texture
- Fractional 1064 HP I (10x10mm), 0.30-0.38J/cm², 5-10hz, 3-5passes
- Fractional 1064 HP II (10x10mm), 0.20-0.25J/cm², 5-10hz, 3-5passes



Before



- Residual Acne, Acne PIH
- Pre-treatment care: deep cleansing with purify peel # 2 sessions in a month
- 1ST : 1064nm, Fractional HP, I 0.30J/cm², II 0.18J/cm², III 0.13J/cm², 7-10hz
- 2ND : 1064nm, Fractional HP, I 0.40J/cm², II 0.30J/cm², III 0.14J/cm², 7-10hz
Glycolic acid 50%(before laser treatment)



Before



After 2Tx

- Pico toning , skin whitening
- 1064 nm, zoom hp, 7mm 10hz, 0.8 J/cm² · 3passes, 2000shots
- 532 nm, S20 hp, 5hz, 0.04 J/cm² , 500 shots




- Wrinkles
- 1064 nm Dia FX HP, 3 sessions per 8 weeks
- Level II, 10x10 mm, 3 hz, 0.24 J/cm²





Periorbital treatment


- Periorbital treatment
- 1st, 2nd TX 1064nm, Dia FX S hp, 5x5mm
- 0.8-1.0 J/cm², 5-10hz, per level I, II, III
- 3rd, 4TH TX, 1064nm, Dia FX S hp, 5x5mm
- 1.0-1.2 J/cm², 7-10hz, per level I, II

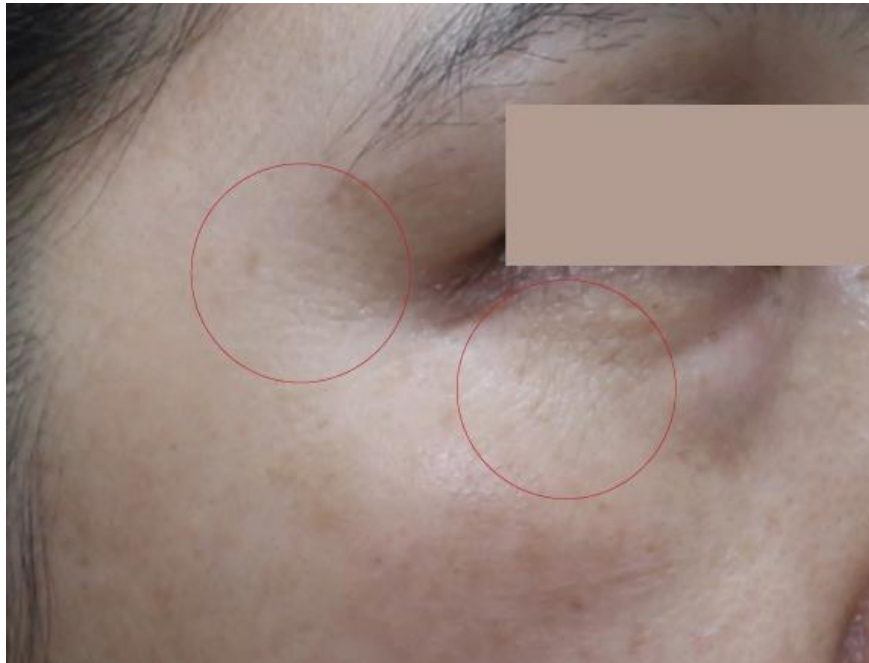

- Periorbital treatment
- 1st, 2nd TX 1064nm, Dia FX S hp, 5x5mm
- 0.8-1.0 J/cm², 5-10hz, per level I, II, III
- 3rd, 4TH TX, 1064nm, Dia FX S hp, 5x5mm
- 1.0-1.2 J/cm², 7-10hz, per level I, II

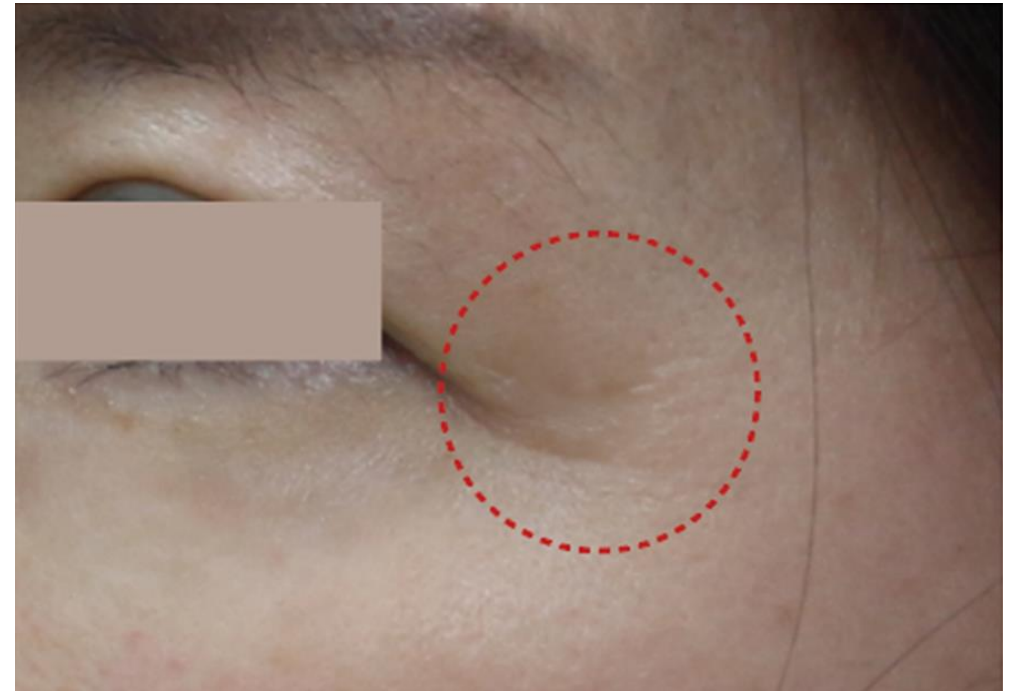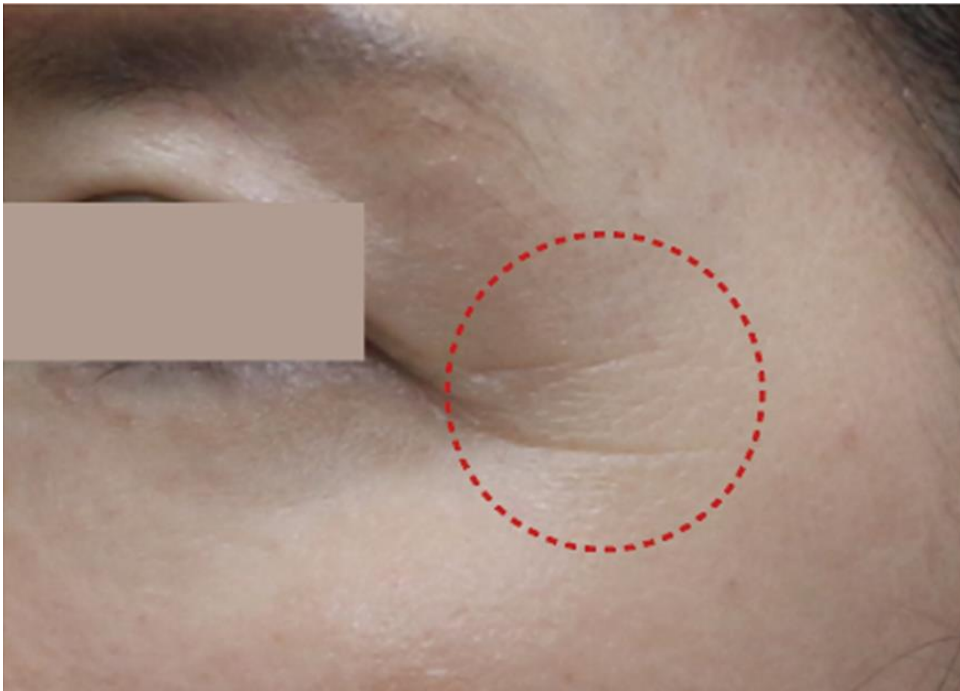

- Periorbital treatment, Pores butterfly zone
- 1064nm, Dia FX S hp, 5x5mm
- 0.8-0.9 J/cm², 10hz, level I, II

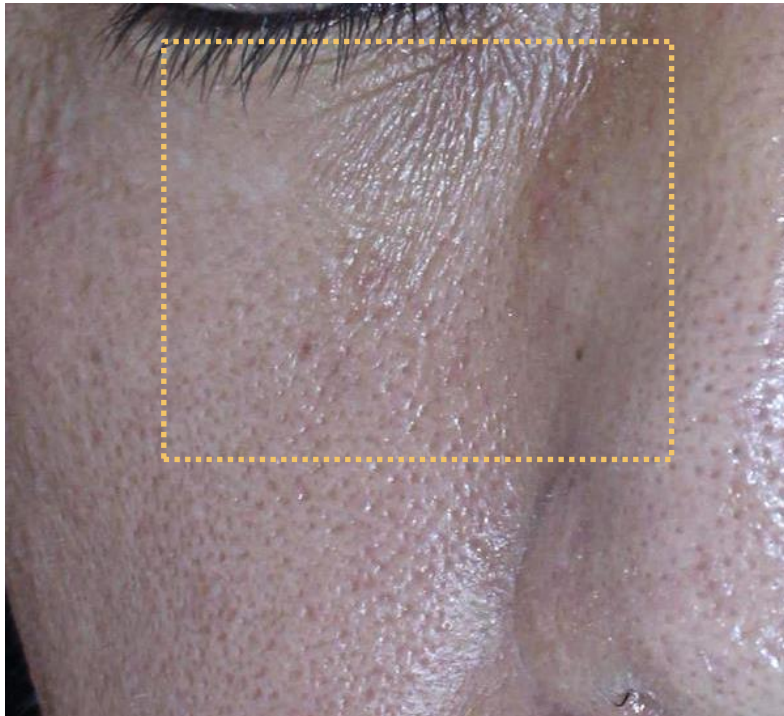


Before




After 3Tx

- Periorbital treatment, Pores butterfly zone
- 1064nm, Dia FX S hp, 5x5mm
- 0.8-0.9 J/cm², 10hz, level I, II

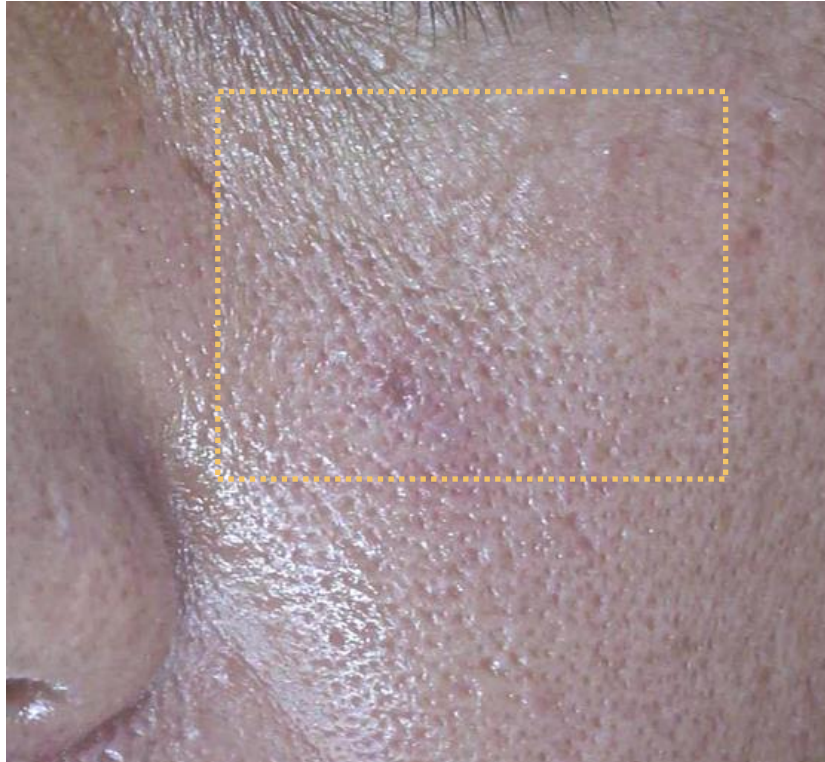

Before

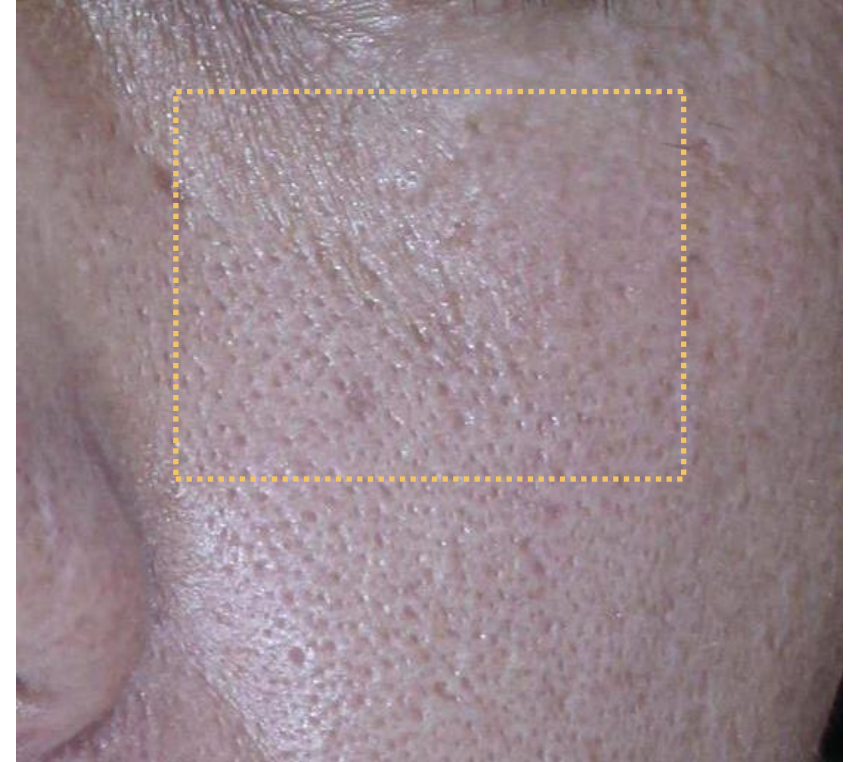

After 3Tx

Pigmentation treatment

- 532nm, Zoom HP, 5mm, 0.6J/cm², 1Hz, 1sessions



Before



After 1Tx

- 532nm, Zoom HP, 5mm, 0.4J/cm², 1Hz, 1sessions




Before



After 1Tx

- Lentigines, Full Pico toning
- Zoom HP, 532nm, 4mm, 0.2-0.3J/cm², 1Hz, 2sessions
- Collimator HP, 1064nm, 10mm, 0.5-0.6J/cm², 10Hz, 8sessions

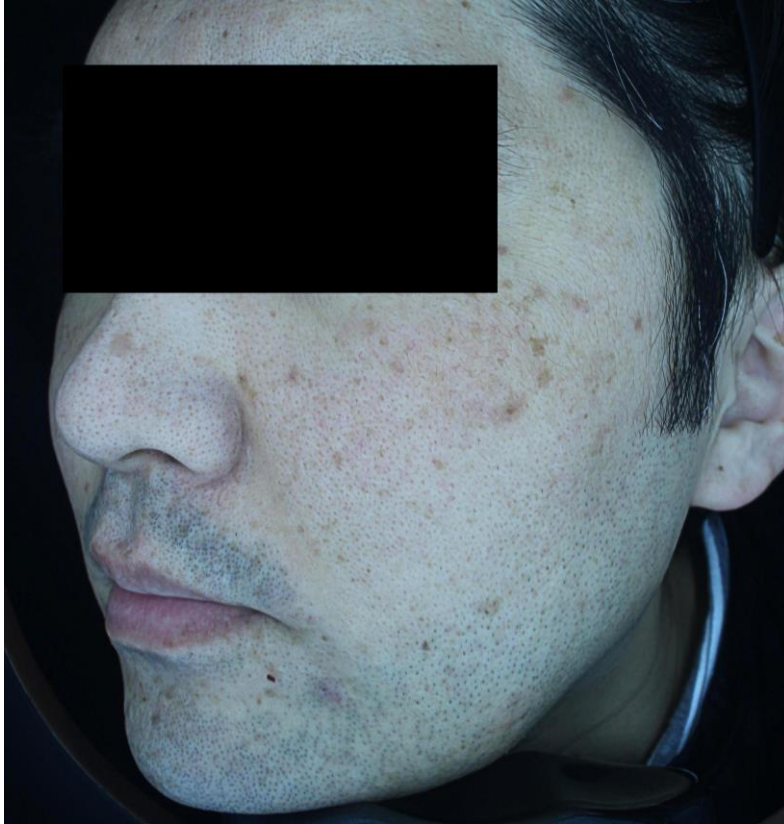


Before




After

- Lentigines, Full Pico toning
- Zoom HP, 532nm, 4mm, 0.2-0.3J/cm², 1Hz, 2sessions
- Collimator HP, 1064nm, 10mm, 0.5-0.6J/cm², 10Hz, 8sessions

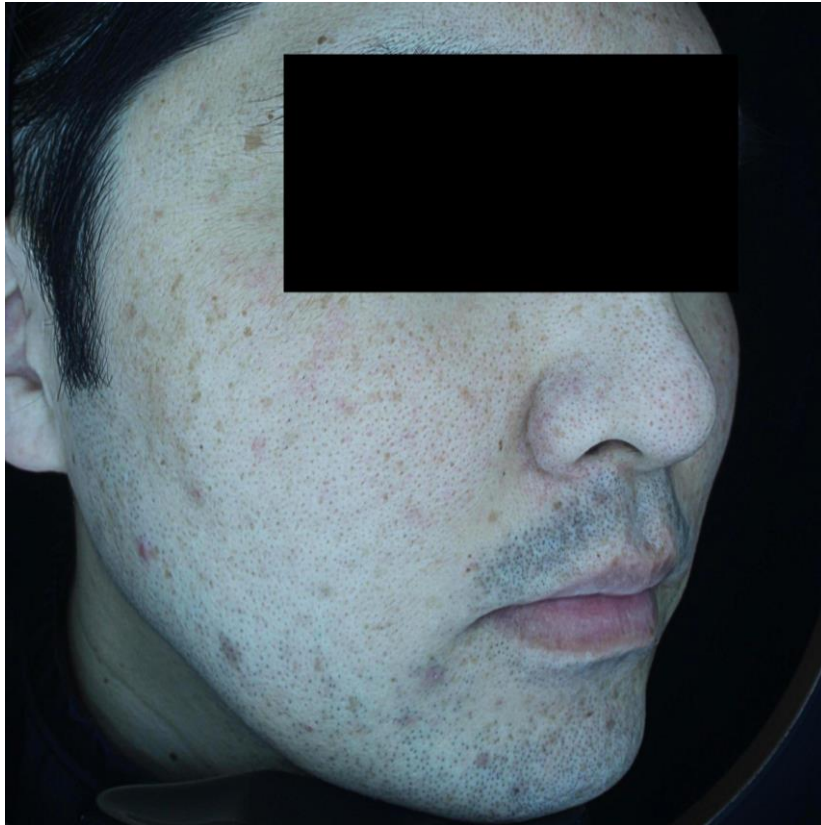

Before

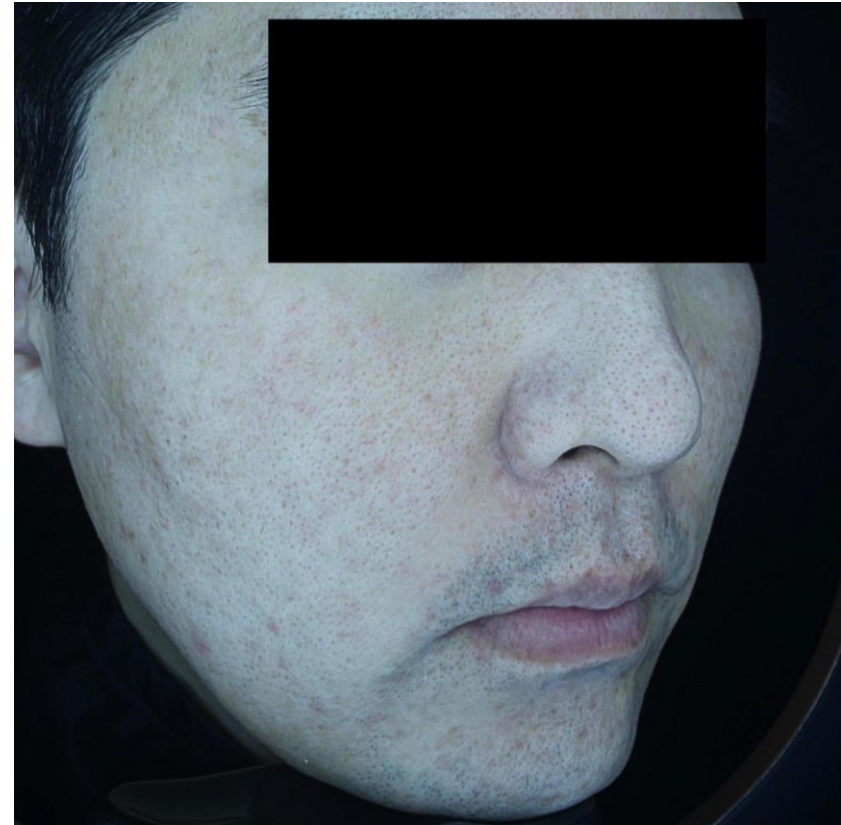

After

- Melasma
- Collimator HP, 1064nm, 10mm, 0.25-0.6J/cm², 10Hz, 3passes



Before



After 6TX

- Lentigo
- Zoom HP, 532 P, 4mm, 0.2J/cm², 1-3hz, 1pass



Before



After 1Tx

- Seborrheic Keratosis
- Zoom HP, 532 P, 4mm, 0.2J/cm², 1-3hz, 1pass



Before



After 1Tx

- Solar Lentigo
- Zoom HP, 532nm with attenuator, 4mm, 1.3J/cm², 1-3hz, 1pass



Before



After 1Tx

- Solar Lentigo
- Zoom HP, 532nm with attenuator, 4mm, 1.3J/cm², 1-3hz, 1pass



Before



After 1Tx

- Solar lentigo
- Zoom HP, 532nm with attenuator, 4mm, 1.3J/cm², 1-3hz, 1pass



Before

After HQ Tx

After PicoLO Tx

- Solar Lentigo
- Zoom HP, 532nm with attenuator, 4mm, 1.3J/cm², 1-3hz, 1pass




Before




Q-switched 532nm
(0.4J/cm²)

After

- Lentigo
- Zoom HP, 532nm with attenuator, 4mm, 1.4J/cm², 1-3hz, 1pass

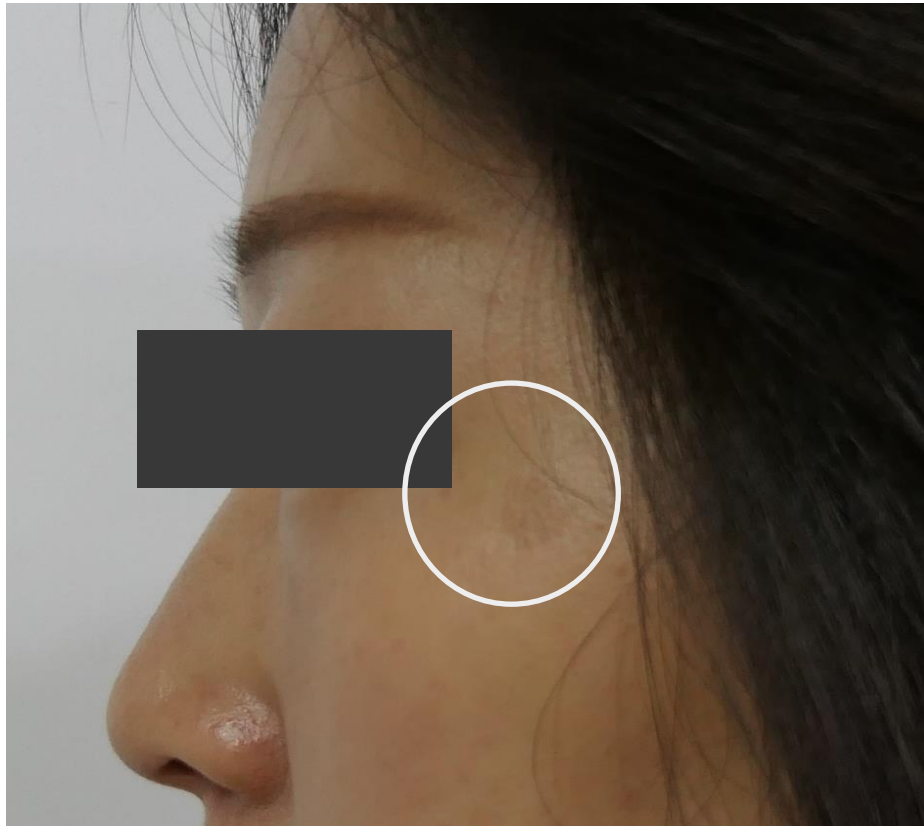

Before

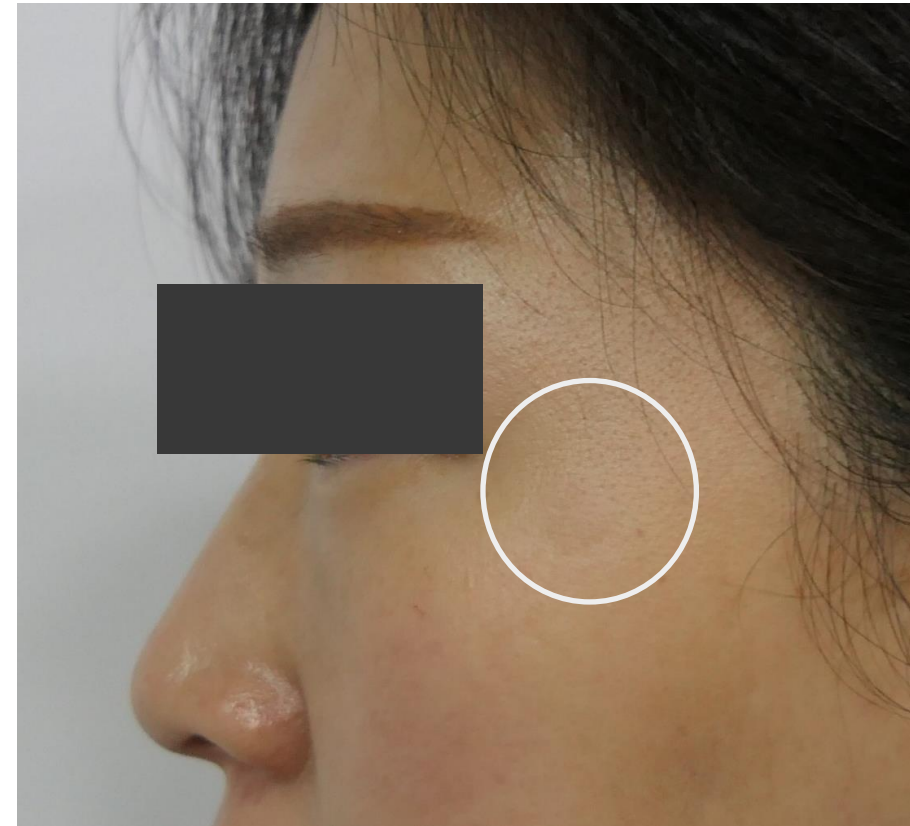

After 1Tx

- Lentigo
- Zoom HP, 532nm with attenuator, 4mm, 1.0J/cm², 1-3hz, 1pass



Before



After 1Tx

- Lentigo
- Zoom HP, 532nm, 4-5mm, 0.5-0.6J/cm², 1-3hz, 1pass




Before



After 1Tx

- Lentigo
- Zoom HP, 532nm, 4-5mm, 0.5-0.6J/cm², 1-3hz, 1pass

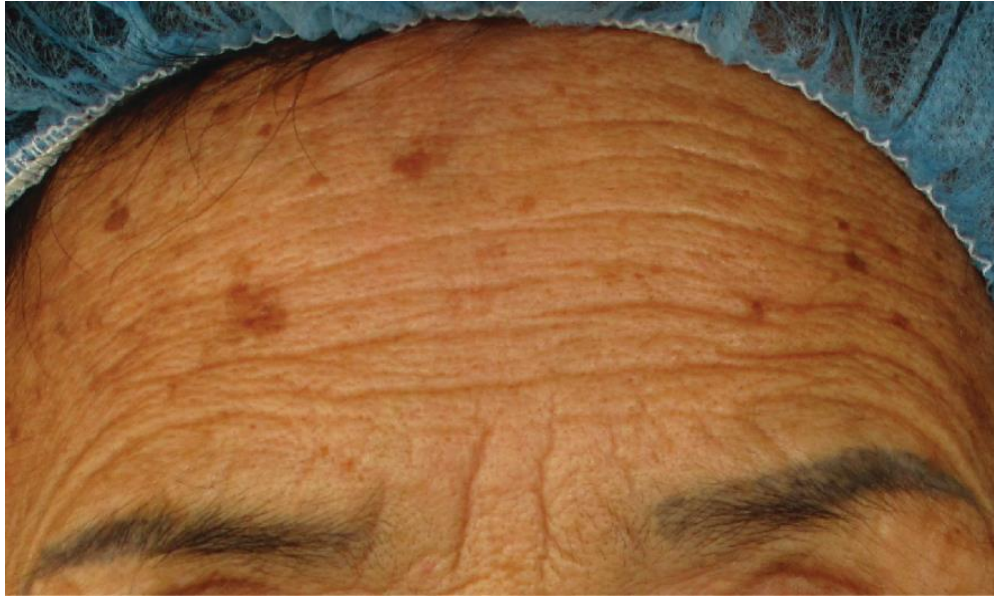

Before



After 1Tx

- Freckles
- 532 zoom 70mJ Zoom HP, 532nm with attenuator, 4mm, 1.0J/cm², 1-3hz, 1pass



Before



After 1Tx

- Nevus of Ota
- Zoom HP, 1064nm, 7mm, 1.0-1.3J/cm², 8-10hz, 3-5passes



Before



After 1Tx



After 2Tx

- Nevus of Ota
- Zoom HP, 1064nm, 7mm, 1.0-1.3J/cm², 8-10hz, 3-5passes



Before



After 5Tx

- Nevus of Ota
- Zoom HP, 4mm, 1064nm, 2.0-4.0J/cm², 1-5hz, 3-5passes




Before



After 2Tx

- Nevus of ota
- Zoom HP, 4mm, 1064nm, 2.0-4.0J/cm², 1-5hz, 3-5passes

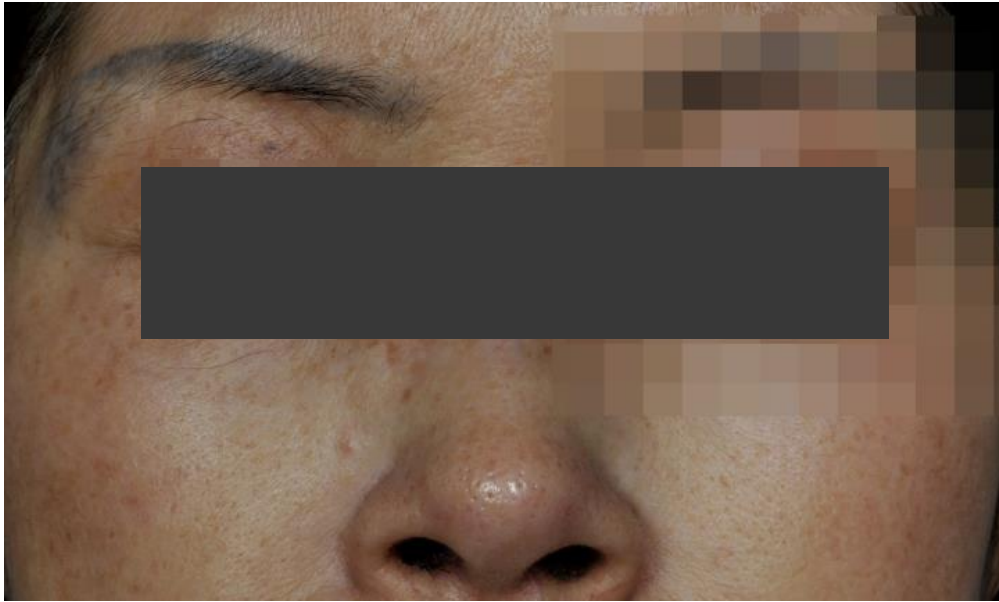

Before



After 2Tx

- ABNOM
- Zoom HP, 1064nm, 4mm, 1.5-2.0J/cm², 5-10hz, 3-5passes



Before



After 5Tx

- FULL TONING
- S20 HP (20x20mm), 1064nm, 0.07J/cm², 8-10hz, 3-5passes



Before



After 7Tx

- Full toning, Melasma, ABNOM, PIH
- S20 HP (20x20mm), 1064nm, 0.07J/cm², 8-10hz, 3-5passes – melasma
- Fractional HP I (10x10mm), 1064nm, 0.30J/ cm², 8-10hz, 3-5passes – PIH



Before




After 7Tx

- Melasma : Zoom HP, 1064nm, 7mm, 0.4J/cm², 8-10hz, 3-5passes
- Ota Nevus, Hori's Nevus : Zoom HP, 1064nm, 3-4mm, 2.0-3.0J/cm², 1-5hz, 3-5passes
- Combination Therapy with Helios III



Before

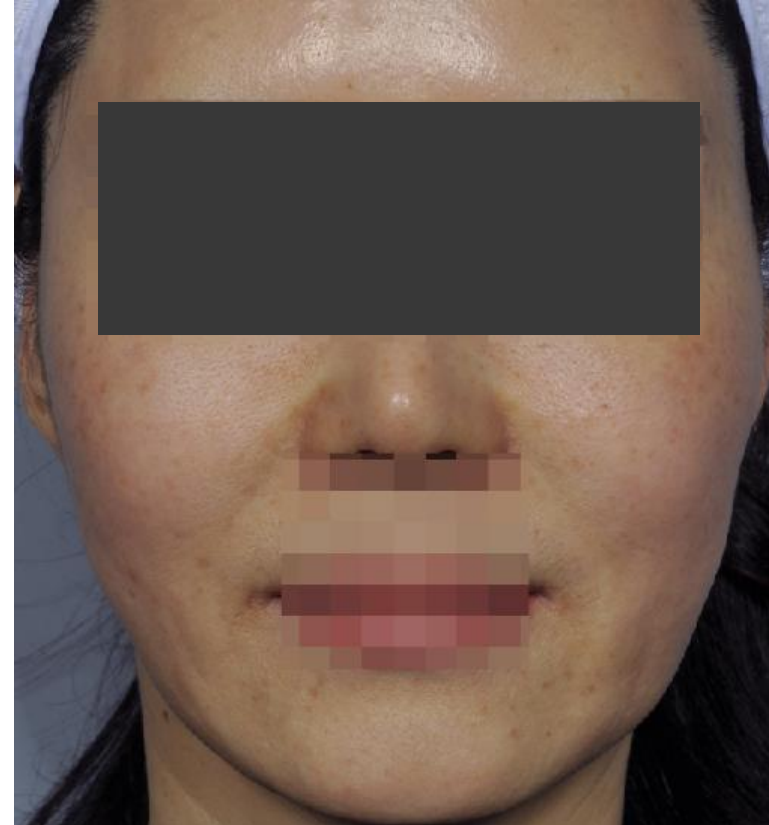

After 2Tx

- Melasma : Zoom HP, 1064nm, 7mm, 0.4J/cm², 8-10hz, 3-5passes
- Ota Nevus, Hori's Nevus : Zoom HP, 1064nm, 3-4mm, 2.0-3.0J/cm², 1-5hz, 3-5passes
- Combination Therapy with Helios III



Before



After 2Tx



- Melanochia
- Fractional 1064 HP (10x10mm), 0.2-0.3J/cm², 5-10hz, 3-5passes

Before



After 5Tx



Before



After 5Tx

- Melanochia
- Fractional 1064 HP (10x10mm), 0.2-0.3J/cm², 5-10hz, 3-5passes

- Solar lentigo
- 532 nm Zoom HP, 1 sessions
- 5mm, 1hz, 0.4 J/cm²



Before



After 1 week



After 4 weeks

- Solar lentigo
- 532 nm Zoom HP, 1 sessions
- 5mm, 1hz, 0.5 J/cm²



Before



Right after treatment



After 4 weeks

- Solar lentigo
- 532 nm Zoom HP, 1 sessions
- 5mm, 1hz, 0.5 J/cm²



Before



After 4 weeks

- Café au lait spot
- 532 nm Zoom HP, 1 sessions
- 4mm, 1hz, 0.6 J/cm²
- Several treatment sessions are required



Before



After falling off crust



After 4 weeks

- Café au lait spot
- 532 nm Zoom HP, 1 sessions
- 4mm, 1hz, 0.6 J/cm²
- Several treatment sessions are required



Before



After 1 week

- Solar lentigo
- 532 nm Zoom HP, 1 sessions
- 5mm, 1hz, 0.4 J/cm²



Before



After 1 week

- Solar lentigo
- 532 nm , Zoom HP with attenuator
- 4mm, 1.3J/cm², 1hz, 1pass



Before




After 4 weeks

- Solar lentigo
- 532 nm , Zoom HP with attenuator
- 4mm, 1.4J/cm², 1hz, 1pass




Before

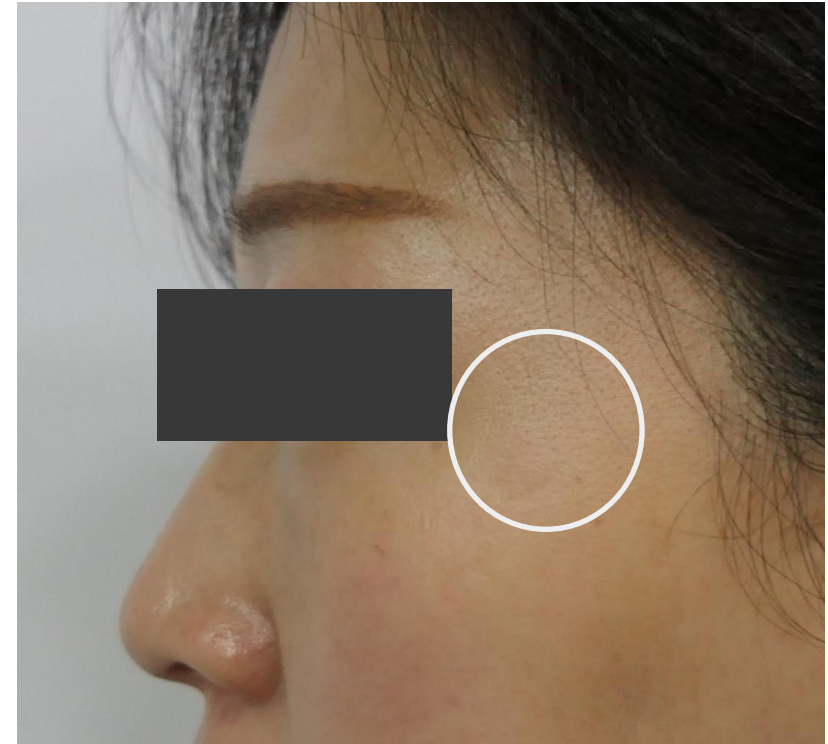


After 4 weeks

- Freckle
- 532 nm , Zoom HP with attenuator
- 3mm, 0.9J/cm², 1hz, 1pass
- Full Pico toning
- 1064nm, Zoom HP, 7mm, 10hz, 1500shots

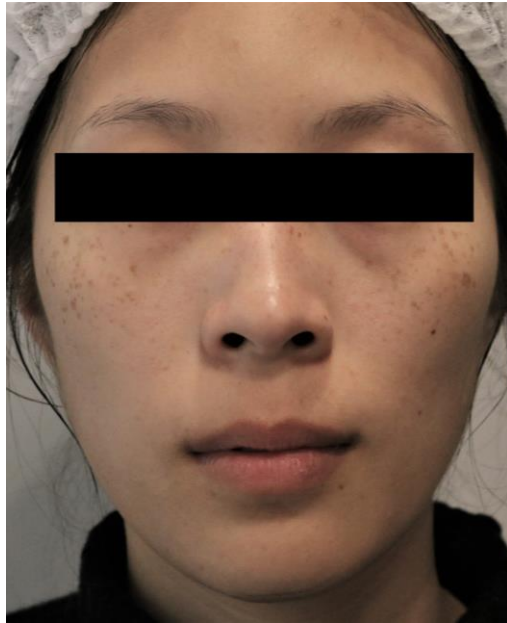

Before

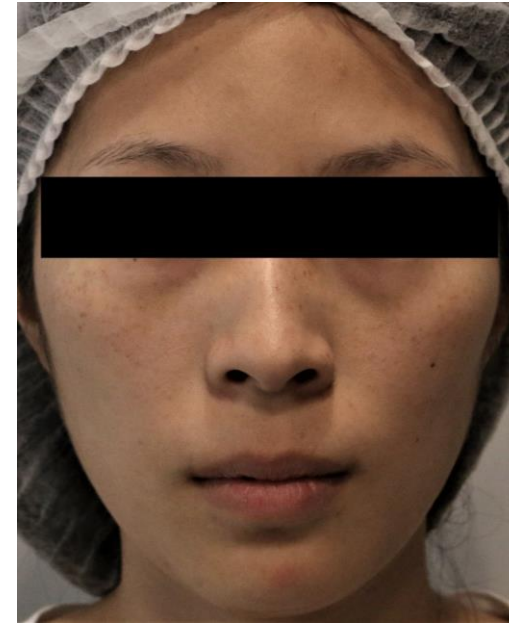


After 4 weeks

- Melasma
- 1064 nm + 532nm S20 HP, 5 sessions
- 1064 nm _ 20x20mm, 10hz, 0.06 – 0.08 J/cm²
- 532 nm _ 20x20 mm, 10hz, 0.02 J/cm²
- Hydroquinone 2%



BEFORE

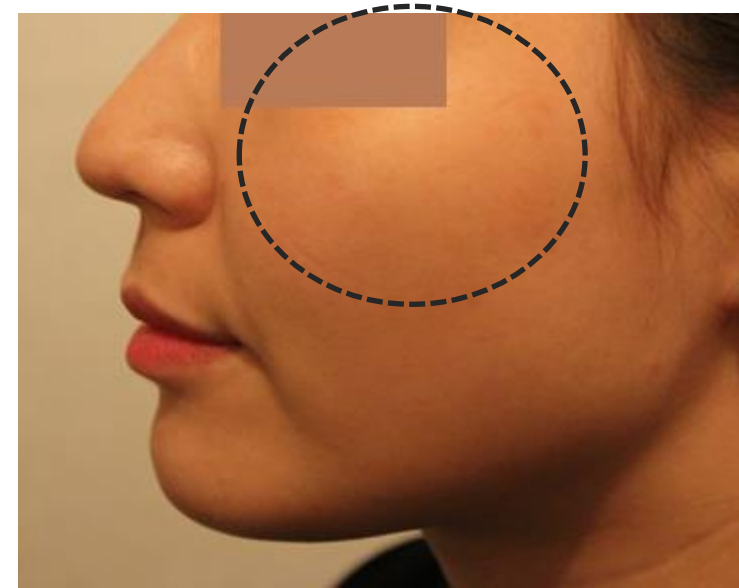

After 5 sessions

- Peutz - Jegher syndrome
- 1064 nm , 3 sessions
- Dia fx, 10x10 mm, 10hz, 0.50 J/cm²



Before



After 3 TX

- Peutz - Jegher syndrome
- 1064 nm , 3 sessions
- Dia fx, 10x10 mm, 10hz, 0.50 J/cm²



Before



After 3 TX

532 Rejuvenation treatment

- Skin rejuvenation, skin whitening, blemishes
- 532 nm, P Mode, 532 Dia FX HP, 2 sessions
- 7x7 mm, 10hz, Attenuator, 0.07 J/cm²



Before



After 2 TX

- Skin tone, skin brightening, rejuvenation
- 532 nm, P Mode, 532 Dia FX HP, 4 sessions
- 7x7 mm, 7-10hz, Attenuator, 0.06-0.1J/cm²



Before



After 4 TX



Before



After 4 TX

- Skin rejuvenation, skin whitening, blemishes
- 532 nm, P Mode, 532 Dia FX HP, 2 sessions
- 7x7 mm, 10hz, Attenuator, 0.08 J/cm²



Before



After 2 TX

- Skin rejuvenation, Pores and Scars
- 532 nm Dia FX HP, 3 sessions
- 7x7 mm, 5-10hz, Attenuator, 0.32-0.44 J/cm²



Before



After 2 TX

- Skin rejuvenation, Pores and Scars
- 532 nm, 1064 nm Dia FX HP, 3 sessions
- 7x7 mm, 10hz, Attenuator, 0.22J/cm²
- 10x10 mm 3hz, 0.15J/cm²



Before



After 2 TX



After 3 TX



Amico Yasna Pars (Pr.J.S.Co)
THANK YOU