

A LEADER is a DEALER in hope.



Inside this issue:

Tecnis Family: 2
Monofocal One Piece
IOLs

Abbott Medical 2
Optics:
Sensar 1-Piece

Huvitz: 3
HNT-7000
Non-Contact
Tonometer

Geuder: 3
megaTRON S4HPS

Abbott Medical 4
Optics:
WhiteStar Signature

Mini Case 4

Vince Lombardi says: "Leaders are not born, they are made and they are made just like anything else, through hard work and that is a price we will have to pay to achieve that goal or any goal."

Dear Valued Partner

There are myriad rules and laws in this world or better say in this cosmos that has been written nowhere. You can just find them by your heart or even by your intuition.

You know why? Because they have been coveted.

One of them is law of attraction. Can you imagine that the energy of a positive

thought has a power 10 times greater than a negative one? It means what you really and positively looking for, you will have it finally. Long story short, thoughts will change to reality. This is the story of **Amico Yasna Pars** staff. We always try to think about good things to make them happen. This is what is common between business and psyche in our plans. Every loyal customer can increase our reputation by 10 times. This is what we eager to see and hear. Let us be your loyal partner as always.

www.amicoyasnapars.com

On behalf of Amico Yasna Pars team

TECNIS® FAMILY OF IOLs

This is Optical Synergy

Abbott
Medical Optics

MONOFOCAL 1-Piece IOLs

SHARPEST VISION

- 70% of patients had best corrected binocular distance visual acuities (BCVA) 20/16 or better with **Tecnis**.
- 96% of patients had best corrected binocular distance visual acuities (BCVA) 20/20 or better with **Tecnis**.
- **TECNIS®** IOLs provide 14–35% improved image contrast from chromatic aberration reduction and they also correct the Spherical Aberration to essentially zero.

ENHANCED FUNCTIONALITY

Improved patient safety under low-visibility conditions

There is likely to be a meaningful safety benefit for elderly drivers with

TECNIS® IOLs and to the drivers and pedestrians with whom they share the road.

Safety that extends beyond driving: the **TECNIS®** IOL improves functional vision, which may improve patient safety for other life situations under low-visibility conditions.

Improved low-light optical performance

TECNIS® 1-Piece Monofocal IOL demonstrates 35% improvement in image contrast with large pupil (5mm).

LONG-TERM SUSTAINABILITY

Unlike another leading IOL, **TECNIS®** IOL material is not


associated with glistenings.

ProTEC 360° barrier edge for reduced LEC migration.



Extracted from:
<http://www.tecnisiol.com>

SENSAR 1-Piece: AMO's very high quality hydrophobic spherical IOL.

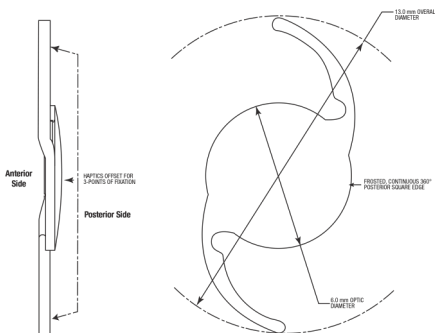

5. Tri-Fix 3-Point Fixation.

6. 360° Posterior Square Edge.

- But unlike **Tecnis** which is Aspheric, The **Sensar** is Spherical IOL which justify its lower price.

Extracted from:
<http://www.amo-distributor-service.com>

Indeed the **Sensar 1-piece** is same as **Tecnis 1-piece** with the same platform which means:

1. Limited LEC Migration.
2. Not Associated with Glistenings.
3. Transmission of Healthy Blue Light.
4. Reduced Chromatic Aberration because of higher Abbe number.

DISCRIPTION

OPTIC CHARACTERISTICS	
Powers	+6.0 D to +30.0 D in 0.5 diopter increments
Diameters	6.0 mm
Material	UV-blocking hydrophobic acrylic
Refractive Index	1.47
Edge Design	ProTEC frosted, continuous 360° posterior square edge
ULTRASOUND BIOMETRY	
A-Constant	118.4
Theoretical AC Depth	5.2 mm



Huvitz



HNT-7000 Non-Contact Tonometer

Huvitz HNT-7000 measures patient's IOP with a simple touch of a button.

Full auto left/right side change, alignment, focusing and puffing is done automatically with incredible ease.

This full auto functionality and corneal thickness data compensation functions provides you with more precise IOP values for diagnosis.

Improved measurement value, better time management and patient comfort are just a few additional advantages of the finest non-contact tonometer from **Huvitz**.

Full Auto Tracking:

Only a button needs to be pressed to align, focus and puff the both, left and right eyes.

High-resolution 5.7" Touch LCD/ User Friendly Interface:

Touch LCD screen with 0 to 90

degree tilting function provides you easy and effective operation for any condition given.

Smart Puffing Control:

Soft puffing maximizes patient comfort. The puffing is automatically regulated after the measurement if needed.

Corneal Thickness Compensation:

For more precise tonometry, corneal thickness values of patient can be inserted to get compensated and accurate data.

With an integrated printer and motorized chin rest, the **Huvitz HNT-7000** combines a sleek and modern design with practicality and patient comfort.

Extracted from:
<http://www.huvitz.com>

Geuder; Be Part of It SOPHISTICATED TECHNOLOGY FOR MAXIMUM PERFORMANCE

The lead is taken by whoever delivers the best performance – and starts from a promising position. Make sure to have the Pole Position. With the sophisticated technology of the new **megaTRON S4HPS** you will achieve a maximum of safety, efficiency, cost effectiveness and flexibility. The revolutionary hybrid pump system is only one of the numerous features which make the new **megaTRON S4HPS** one of the most powerful device systems for the anterior and posterior segment. In combination with

the UNO Colorline MACH2vitreous cutters from now on cut rates of up to 12,000 cpm are possible. The optional endolaser **endoTRON532** and the xenon light-source **xenoTRON III** complete the **megaTRON S4HPS** and make it an all-rounder in vitrectomy. Of course we provide you with a large variety of accessories for the anterior and posterior segment – perfectly adjusted to the complete **megaTRON S4HPS**.

Extracted from:
<http://www.geuder.de/>

Geuder®
Precision made in Germany





Advanced Medical Optics (AMO) positions its new **WhiteStar Signature** with Fusion Fluidics as offering flexible new capabilities to

WhiteStar Signature: Next-generation phaco increases chamber stability with precise fluidics control

minimise surgical trauma by reducing phaco energy and post-occlusion surge. This is done by integrating advanced fluidics management that allows higher vacuum and flow with greater chamber stability, and enhanced micropulse technology that increases emulsification efficiency. The system is designed to support both existing and anticipated coaxial and bimanual procedures using a variety of current and emerging surgical techniques, particularly **micro-incision** approaches. **Signature** is designed for efficiency in the OR, featuring a compact footprint and access to pre-programmed fluidic and **phaco** power settings through a

touchscreen interface. Snap-in tubing cassettes enable quick turnaround between cases with reduced chances for set-up errors. The heart of the system is a dual pump that provides both venturi-type vacuum and peristaltic flow controlled by Fusion Fluidics management. Fusion Fluidics incorporates what **AMO** calls advanced chamber stabilisation environment, to reduce post-occlusion surges. "It's like ABS brakes, you don't have to think about it," said David F Chang, MD, clinical professor, University of California San Francisco, US, who consults with **AMO**.

Extracted from:
<http://www.es CRS.org/>

Mini Case:

A 22-year-old man had perceived a moving black spot in the visual field of his right eye and no other symptoms for 1 year. An ophthalmologic examination revealed visual acuity without correction of 20/20 in both eyes, and slitlamp examination findings of the right eye identified a unique lesion that was oval in shape, measured 2 mm in thickness by 2 mm in diameter, was translucent and brown, and changed its location with head posture.

During gonioscopy, we observed a nonpigmented wide angle without other findings that suggested adhesions between the lesion and the structures of the cameralar angle. The patient had a clear

cornea without pathologic findings, and the intraocular pressure with Goldman tonometer was 10 mm Hg in both eyes. Findings from a dilated fundus examination were unremarkable

A Moving Black Spot in My Vision: What would you do next?

- A computerized tomographic scan
- Anterior segment tomography with Pentacam
- Begin systemic treatment with albendazole
- Surgical removal of the cyst and pathological analysis

You can also refer to <http://www.archophth.jamanetwork.com>

