



TECNIS Synergy™ IOL

Johnson & Johnson
Amico Yasna Pars (Pr.J.S.Co)

Get the best of both worlds.
Focus far through near;
clarity day and night.

Widest range with best near

Superior image contrast

Excellent patient outcomes

An innovation that combines the best of TECNIS® Multifocal and TECNIS Symphony™ technologies

TECNIS Synergy™ IOL combines diffractive technologies

TECNIS® Multifocal IOL technology


- Good near VA⁷
- No compromises on distance vision⁷
- Pupil-independent design⁷

TECNIS Symphony™ ERV/EDOF IOL technology

- Continuous vision across the range⁸
- **ChromAlign™** Technology
- Pupil-independent design^{8,9}

The proprietary achromatic technology you know and trust from the TECNIS Symphony™ IOL is now available in the TECNIS Synergy™ IOL


ChromAlign™ Technology corrects corneal chromatic aberration for enhanced image contrast across the range of vision^{8,10}

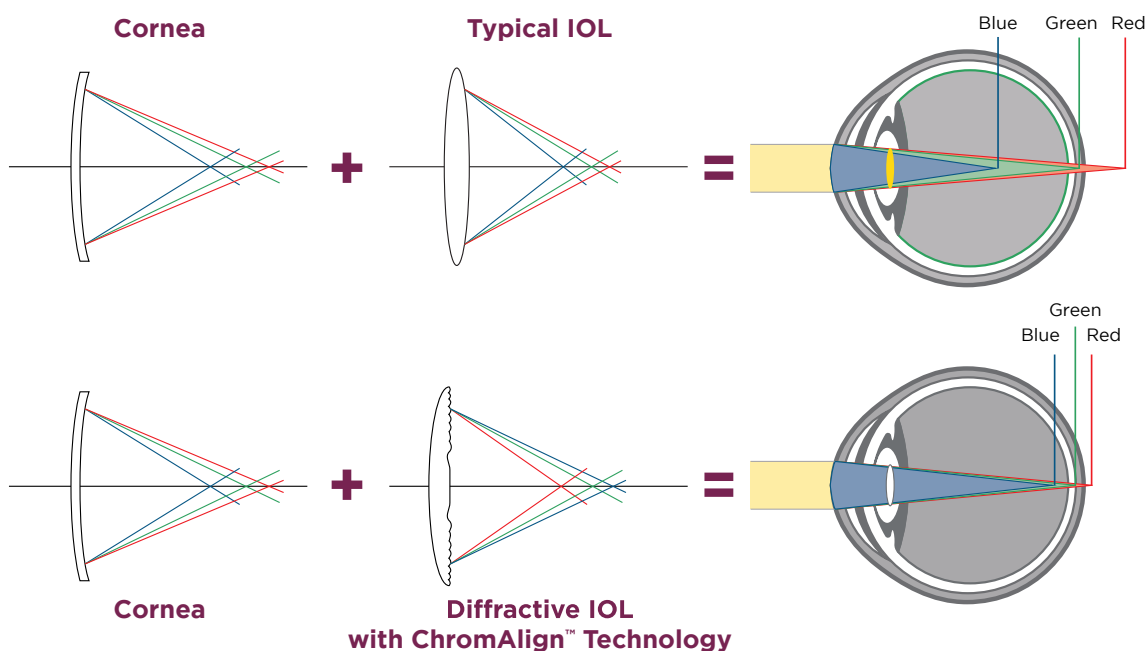

TECNIS
Synergy™ IOL

Continuous vision that empowers®
patients to see at their ideal distances

TECNIS Synergy™ IOL: Unparalleled “in the zone” vision that outperforms trifocal IOLs, including outstanding near visual performance at 33 cm.

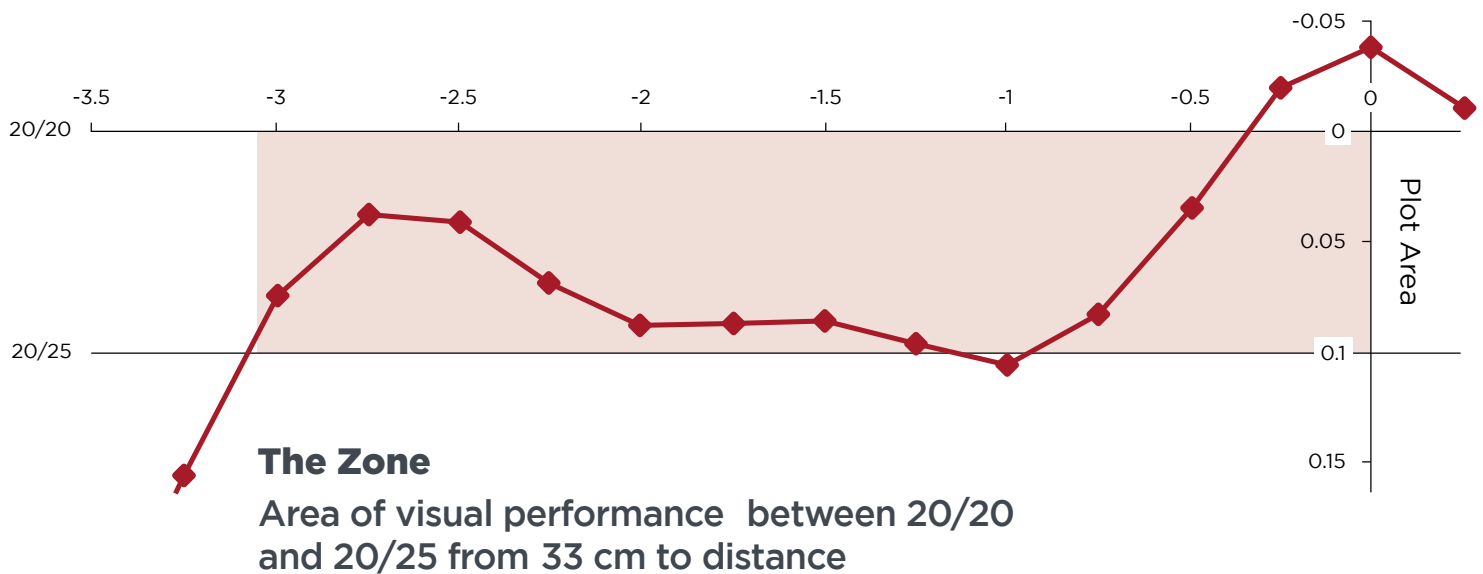

Violet-light filtration reduces visual symptoms such as halos

With violet-light filtration, ease of both daytime and nighttime driving showed statistical improvement



With violet-light filtration, halo intensity was reduced in both iPhone screens (photopic -13%) and xenon car headlights (scotopic -29%)



OPTICAL CHARACTERISTICS¹

Powers:	+5.0 D to +34.0 D in 0.5 diopter increments
Diameter:	6.0 mm
Center Thickness:	0.7 mm (20.0 D)
Shape:	Biconvex, wavefront-designed anterior aspheric surface and achromatic technology to correct chromatic aberration for enhanced contrast sensitivity.
Material:	Soft, foldable hydrophobic acrylic with UV and violet light absorber
Refractive Index:	1.47 at 35° C
Edge Design:	ProTEC frosted, continuous 360° posterior square edge
Achromatic Technology:	Proprietary technology for chromatic aberration correction

BIOMETRY¹

CONTACT ULTRASOUND¹

OPTICAL²

A-constant:	118.8	119.3
AC Depth:	5.4 mm	5.7 mm
Surgeon Factor: ²	1.68 mm	1.96 mm

HAPTIC CHARACTERISTICS¹

Overall Diameter:	13.0 mm
Thickness:	0.46 mm
Style:	C, TRI-FIX haptics offset from optic; 1-piece lens
Material:	Soft, Foldable, UV-absorbing and violet-light filtering hydrophobic acrylic

Preloaded **TECNIS Simplicity**® Delivery System



Please scan the
QR code to see
instruction for use

BENEFITS:



Easy 3-step process -
hydrate, advance, deliver



Advanced design for a
smooth and controlled
IOL delivery



TECNIS Simplicity™ — The Next-Generation
TECNIS® IOL Disposable Delivery System

Simplifying delivery of the IOL you prefer.

Johnson & Johnson
Amico Yasna Pars (Pr.J.S.Co)

Tel: +9821-71137000

Fax: +9821-71137300

2nd Floor, No. 1698, Shariati Ave., Tehran, Iran

Postal Code: 1914744755

www.aypmedical.com

info@aypmedical.com

